DAAAM INTERNATIONAL VIENNA

www.daaam.com



ANNALS OF DAAAM FOR 2011 & PROCEEDINGS

OF THE 22ND INTERNATIONAL DAAAM SYMPOSIUM
"INTELLIGENT MANUFACTURING & AUTOMATION:
POWER OF KNOWLEDGE AND CREATIVITY"
23-26TH NOVEMBER 2011, VIENNA, AUSTRIA

ORGANIZED BY:

DAAAM INTERNATIONAL VIENNA
INTERNATIONAL ACADEMY OF ENGINEERING
VIENNA UNIVERSITY OF TECHNOLOGY, UNIVERSITY OF APPLIED SCIENCES
TECHNIKUM VIENNA AND AUSTRIAN SOCIETY OF ENGINEERS AND ARCHITECTS - ÖIAV 1848

Under the auspices of: The Danube Rectors' Conference & Rectors' Honor Committee of DAAAM International

EDITOR: B.[RANKO] KATALINIC







Verba volant, littera scripta manet. [In English: Words fly away, written letters remain.]





Jubilee Gold Medal of DAAAM International

Who is in the field, and not in the DAAAM he/she does not exist. (Jyri Papstel 1996)

EDITOR'S NOTE This publication was published using the soft copy manuscripts supplied by their authors. Manuscripts are subject of the following review procedures: (a) It is reviewed by editor for general suitability for this publication. (b) If it is judged suitable two reviewers are selected and a double blind review process takes place. (c) Based on the recommendations of the reviewers, the editor then decide whether the particular contribution should be accepted as it is, rev sed or rejected. The editor and the DAAAM International Vienna are not responsible either for the statements made or for the opinion expressed in this publication.

Copyright © 2011 by the DAAAM-2011 editor Abstracting and non-profit use of the material is permitted with credit to the source. Libraries are permitted to photocopy for private use of patrons. Instructors are permitted to photocopy isolated articles for non-commercial classroom use without fee. Author(s) agree that the DAAAM international editor 2011 (B[ranko] Katalinic) receives a full and exclusive, worldwide, perpetual, non-revocable and transferable license under copyright in the submitted work. The licence especially covers, without limitation, all rights to make the article available in the internet or databases or to cite to reference, to cross-reference, to publish, reproduce and distribute the article. However the author(s) are entitled to make the following use of their work: a) to republish the article in any publication of which they are an author or editor after the work has been published by the DAAAM editor; b) to photocopy the article and distribute such photocopies and distribute the electronic versions of the article for personal or professional use, provided that any such copies are not offered for sale, and c) to make oral presentations of the work in any forum. Any of the above mentioned use of the work must explicitly identify prior publication in the *Annals of DAAAM for 2011 & Proceedings of the 22nd International DAAAM Symposium, ISSN 1726-9679, ISBN 978-3-901509-83-4, Editor B[ranko] Katalinic, Published by DAAAM International, Vienna, Austria 2011, including page numbers.*

DAAAM International is a member of CrossRef. Contributions published in DAAAM Scientific Book since 2006 and DAAAM International Journal of Simulation Modelling are cross-referenced and available on Internet. DAAAM International publications are cited and referenced in the: ISI Scientific Proceedings is now the Conference Proceedings Citation Index (CPCI) - an integrated index within Web of Science Thomson Reuters. EBSCO Publishing - bibliographical database. Papers are available at Thomson Gale.

ANNALS OF DAAAM FOR 2011 & PROCEEDINGS OF THE 22ND INTERNATIONAL DAAAM SYMPOSIUM, Vol. 22, No1, Abbreviated key title: Ann. DAAAM proc. int. DAAAM Symp.

ISSN 1726-9679, ISBN 978-3-901509-83-4, Editor B[ranko] Katalinic Published Annually by DAAAM International Vienna, Vienna 2011

Layout & Design: B[ranko] Katalinic / Technical Editors: P[aulina] Cesarec / R[oman] Kettler / I[lya] Kukushkin DAAAM Home Page Admins: V[iktorio] Malisa / M[arkus] Stopper / D[imitrios] Prodromou / Domagoj Kriskovic Publisher: DAAAM International Vienna, TU Wien, Karlsplatz 13/311, A-1040 Vienna, Austria

E-mails: president@daaam.com katalinic@mail.ift.tuwien.ac.at

DAAAM homepage: http://www.daaam.com http://www.daaam.org http://www.daaam.info Graphics used in this Publication: M.C. Escher



Annals of DAAAM for 2011 & Proceedings of the 22nd International DAAAM Symposium, Volume 22, No. 1, ISSN 1726-9679
ISBN 978-3-901509-83-4, Editor B. Katalinic, Published by DAAAM International, Vienna, Austria, EU, 2011

Make Harmony between Technology and Nature, and Your Mind will Fly Free as a Bird

Annals & Proceedings of DAAAM International 2011

APPLICATION OF THE SIMULTANEOUS DSC/DTA/QMS THERMAL ANALYSE TO ALLOY AND NANOPOWDER CHARACTERIZATION

PESINA, Z[bynek] & SOPOUSEK, J[iri]

Abstract: The unique simultaneous thermal analyse is introduced. The individualy adapted STA 409 / skimmer system deals with possibilites of application of thermal analysis methods. Methods of thermal analysis (DTA, DSC, QMS) were used to determine a wide range of thermodynamic quantities (e.g. melting point) of pure metal and alloys. DSC experiments were carried out for both bulk materials and for nanopowders. The paper also indicated the potential of CALPHAD method. Key words: DTA, DSC, QMS

1. INTRODUCTION

Thermal analysis is one of the most important methods of studying the composition and physic-chemical properties of materials. Due to dynamic development of this method there are many manufactures in the market (Setaram, TA Instruments, Netzsch and Linseis). Methods of thermal analysis are useful in studying different samples and monitoring their utility properties. These devices are usually equipped with options to operate different methods of thermal analysis and ordinarily associated with other analytical methods (thermogravimetry TG, mass spectrometry MS, interfered spectroscopy IR, Fourier transform infrared spectroscopy FTIR). The materials examined with the aid of thermal analysis can be very diverse.

Commercially available products are adapted to customer requirements. These devices are particulary suitable for use in basic research. This present contribution introduces one of these exceptional instruments. Its uniqueness lies in a connection of simultaneous thermal analysis with quadrupole mass spectrometer to STA 409 without skimmer system. This adaptation enables application of holder with Knudsen cell. The resulting configuration is capable of detecting extra-low gas pressures in metals and alloys.

2. EXPERIMENTAL APPARATUS

The simultaneous thermal analysis was carried out using individually adapted Netzch STA 409 skimmer apparatus for thermal analysis coupled with quadrupole mass spectrometer (QMS). The device enables to apply following instrumental methods: heat flow differential scanning calorimetry (DSC), differential thermal analysis (DTA) and Knudsen cell combined with mass spectrometry (KC/MS). The temperature range is (25-1400) °C, rate of heating and cooling is within the range of (0.1–25) K/min. The gas (argon, nitrogen, air) flow range is (0-300) ml/min. The results measured can be treated via Protheus software. The different crucibles (Alumina, Pt) covered with protective layers can be used with respect to system under view.

3. DSC

The high sensitivity DSC was confirmed on various systems (e.g. Sn, Ag, Cu, Fe Ni, Al). The optimum of sample is

approximately (10-30) mg. The high sensitivity of DSC measurement at the heating rate of 0.1 K/min is demonstrated by example (lead free alloy) in Fig. 1. Reduction of heat rate allows separating a complex signal, which is formed by overlapping large peaks and shoulders. The system is more close to equilibrium state and the transformation temperatures of the samples can be attained with approximately 0.3 K precision. Also, the high sensitivity enables also measurement of the heat effects of the transformations solid to liquid as well as in the solid state. The heat capacity measurement is also possible.

The better understanding of the experimental DSC signals can be achieved in connection to CALPHAD based method (Boettinger et al., 2006) compare with Fig.2. The CALPHAD approach enables each alloy studied to obtain the predicted temperature in dependence on the enthalpy function (Δ H) and the values of temperature enthalpy derivatives $\partial \Delta H/\partial T$ can consequently be used for comparison with DSC experimental signal.

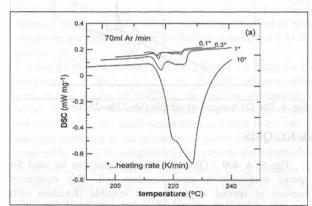


Fig. 1. Experimental DSC heating curves of sample Sn-1.0Ag-0.5Cu-1.0Bi at different heating rates (Palcut et al., 2009)

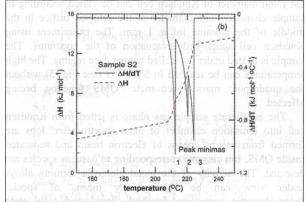


Fig. 2. Calculated molar enthalpies of the Sn-1.0Ag-0.5Cu-1.0Bi and the respective numerical temperature derivative (Sopousek et al., 2009)

The DSC method can be used for bulk alloys as well as for nanopowders, see Fig. 3. The melting point depression, which is the best known property of nanoalloys can be easily detected.

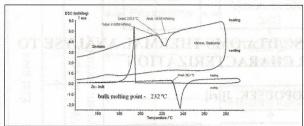


Fig. 3. The DSC signal for bulk and nanoscale tin (Broz et al., 2009)

4. DTA

After a simple holder change, the device Netzch STA 409 offers the possibility of running of the DTA method. For DTA samples both alumina crucibles and silica ampoules (temperature limit of silica being 1100 °C) can be used. The method is suitable for determination of the phase transformations temperature of alloys. However, larger amounts of sample (approximately 0.5 cm³) than in the DSC method are needed. In sample sensitive to oxidation or high vapour pressure, using silica ampoules is unavoidable. Thus, the mass loss of sample is prevented and the observed signal is highly reproducible, see Fig. 4.

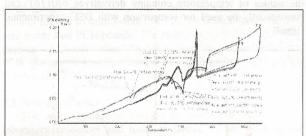


Fig. 4. The DTA signal of sample (Sb-23Sn-25Zn)

5. KC/QMS

The STA 409 / QMS (403/5) apparatus can be used for special measurement using Knudsen cell. This equipment consists of special single alumina crucible (Knudsen cell) covered by lid with small orifice. This crucible is placed on holder top with thermocouple. This arrangement can easily alternate with DSC or DTA holders. The volume (approx. 5 cm³) and diameter of Knudsen cell can be varied according to sample character. The diameter of the effusion orifice in the middle of the alumina lid is 1 mm. The experiment using Knudsen cell starts with evacuation of the apparatus. The sample is heated under controlled temperature regime. The high temperatures can be reached in STA 409 (see Fig. 5) without the quadrupole mass spectrometer QMS function beeing affected.

The equilibrate sample gas phase is jetted from Knudsen cell into ionisation chamber of mass spectrometer. Ions are formed from gas by means of electron beam and separated inside QMS. Ion currents corresponding to gaseous species are detected. The activity of the metal elements forming alloys under view can be obtained by means of special thermodynamic procedure for both liquid and solid state measurements. The example of the activity determination is given in Fig. 6 (Sopousek et al., 2010)

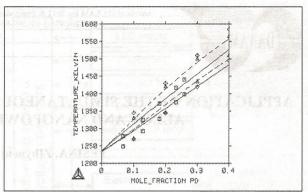


Fig. 5. Liquid and solidus temperatures calculated or predicted by different authors (lines) and selected experimental values (symbols)

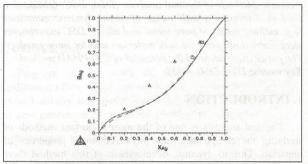


Fig. 6. Activity Ag in liquid AgPd system at 1873 K. Predicted dependences – lines, experimental values - symbols

6. CONCLUSIONS

DSC and DTA methods enable the temperature of phase transformations to be measured with high accuracy. DSC method is more suitable for heat effect evaluations. DSC and QMS measurement using Knudsen cell can be used for special applications as thermodynamic activity evaluations for metal element in alloys.

7. ACKNOWLEDGEMENTS

Financial support from the Czech Science Foundation (grant No. 106/09/H035) and Ministry of Education, Youth and Sports of the Czech Republic (MSM0021622410, CEITEC) is greatefully acknowledged.

8. REFERENCES

Boettinger, W.J. & Kattner, U.R. & Moon, K.W. & Pererezko, J.H. (2006). DTA and DSC Measurement of Alloy Melting and Freezing. Washington, National Institute of Standards and Technology, 2010

Broz, P. & Sopousek, J. & Vrestal, J. (2009). On the stability of Ag, Cu and Sn nanoparticles. In CALPHAD XXXVIII, Prague 17-22.5. 2009, pp. 118-118

Palcut, M. & Sopousek, J. & Trnkova, L. & Szewczykova, B. & Hodulkova, E. & Ozvold, M. & Turna, M. & Janovec, J. (2009). Thermal analysis of selected tin-based lead-free solder alloys. *Metallic materials* 47/1, pp. 43-50

Sopousek, J. & Palcut, M. (2009). The interpretation of the DSC signals of the lead-free solder alloys using the CALPHAD approach. In CALPHAD XXXVIII, Prague 17-22.5. 2009, pp. 167-167

Sopousek, J. & Zemanova, A. & Vrestal, J. & Broz, P. (2010). Experimental determination of phase equilibria and reassessment of Ag-Pd system. *Journal of Alloys and Compouds*. Vol.504., No.2, 2010, pp. 431-434